



Visual Call 5000 System

System Description

The *CORNELL* 5000 Visual Call System is the standard for simplicity in terms of design, installation, operation and maintenance. More than a cost effective, easy-to-use Visual Nurse Call System, the 5000 Series offers customers many choices of custom configuration, thanks to the comprehensive field programming capability built into the software. Serial Data lines control the addressable Patient Stations, Auxiliary Stations, Dome Lights and Annunciator Panels. On-site custom configuration allows selection of the system layout and size, station types, zone layout, and output devices for each addressable station. The terminal may also be used to display the status of to each Patient Station for emergency response and diagnostic purposes.



Benefits

- Call Levels: *Staff Present, Routine, Emergency, Assist, Code Blue, Auxillary*
- 32, 64 or 128 addressable Patient Stations per system
- Duty Stations – 8 for each patient station channel
- Nurse locator function available
- Annunciator Panels – 4 maximum each handling up to 32 stations
- Room mapping, zoning and custom configuration on-site
- Lamp testing of all dome lights and LEDs
- Continuous supervision of all addressable devices
- Membrane style call and reset switches – easy to clean and maintain
- Quick-connect/disconnect field wiring terminations – for easy installation
- Service to eight-conductor corridor cable



Pull Cord Emergency Station





Digital Staff Console

5000 series - continued



Visual Annunciator




Patient Station

Common Model Identification Numbers

Visual Annunciator Panel, 32 Station, 4-gang	HRN7614UA	Call Cord, Single	HAN7608UA
Digital Annunciator Console	HCN 7636/38	Satellite Stations	HWN7602/03/04UA
Pull Cord Emergency Station, 1gang	HWN7605/52/53UA	Economy Satellite Station	HWN7624/25AA
Check-in/Status Station	HWN7726UA	Duty Station 3 LED, 2-gang	HWN7606UA
Pushbutton Emergency Station, 1gang	HWN7602UA	Smart Stations, 1gang	HWN7628UC
Staff Locator Station, 1gang	HWN7603UA	Dome or Zone Light, 1, 2, 3, 4 or 5 lamp	HWN669XUA
Code Blue Station, 1gang	HWN7604UA	Control Unit, 32 Station	HTN7616UA
Patient Stations, 1gang	HWN7600/7630UC	Control Unit, 64 Station	HTN7617UA
Economy Patient Stations, 2 & 3 gang	HWN7668/78UC	Control Unit, 128 Station	HTN7680UA



Pillow Speaker

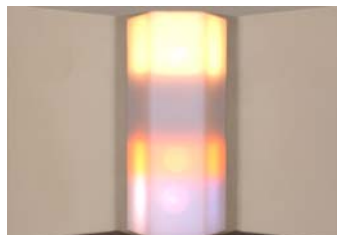

Dome Light



Bedside Station

